

MIS® | LOCKIT

An Affordable, Long-Term, Quality Solution

An Affordable, Long-Term, Quality Solution

The new MIS LOCKiT provides advanced anchoring for implant-retained overdentures. The system was designed to answer the need for a long-term solution, while maintaining a high level of quality, patient comfort and affordability.

The system was designed with customizable levels of retention and features a concave emergence profile, anodized color-coding for easy identification during workflow and is available for MIS internal hex. and conical connections, in all platforms; narrow, standard and wide.



Worry-free reliability

The LOCKiT is coated with anti-wear material and is strengthened above the threads, which enables maximum durability. The plastic caps exhibit better resistance to abrasion over time and have a more predictable and consistent retention level.



Versatility

The extended range of components enables restoration of a non-parallel implant with up to 20 degrees of angulation. This adds up to an extensive 40 degrees of divergence between two implants.



Esthetics

The LOCKiT attachment was designed to reduce the height of the attachment for endosseous implants. The 1.5mm height of the attachment makes this one of the shortest on the market, which provides a convenient and esthetic solution for the patient.



Tissue preservation

The concave emergence profile of the LOCKiT was designed to enable more room for soft tissue growth. Higher volume of soft tissue may potentially lead to bone preservation.

Lockit System Workflow

The following illustrates the recommended workflow for the LOCKiT system placement procedure.



① LOCKiT abutment selection




② LOCKiT installation




③ Silicon rings and metal housing


note : In cases where the concave profile emerges above soft tissue margin, a wider silicon ring (MB-DB235, included in the kit) may be used.

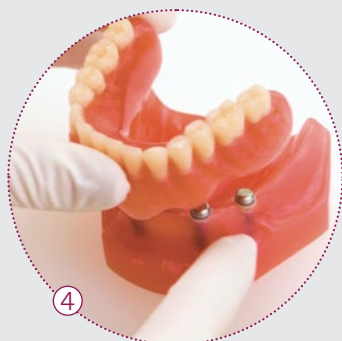


④ Denture preparation

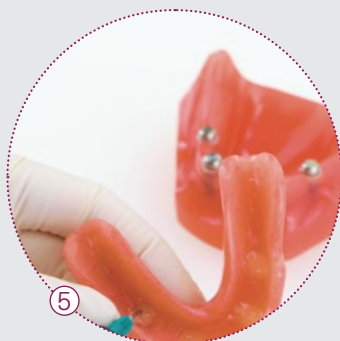


⑤ Denture re-line

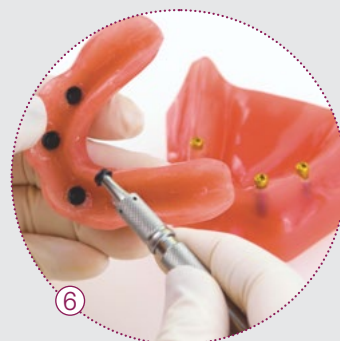

⑥ Male processing cap removal

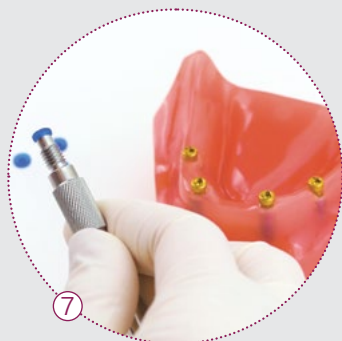

⑦ Male retention cap selection

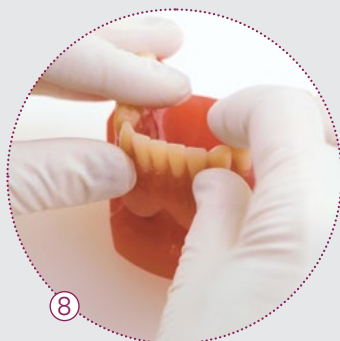




⑧ Try-in and delivery







⑨ Follow-up visits



























Item List


LOCKiT Analog & impression (4-pack)		
MM-LOA40	LOCKiT analog Ø4mm, NP/SP	 \$68
MM-LOA50	LOCKiT analog Ø5mm, WP	 \$68
MM-LOC11	LOCKiT impression coping	 \$68

LOCKiT plastic caps (4-pack)		
MM-CLE10	LOCKiT male retention cap (clear)	 \$20
MM-PIN10	LOCKiT light retention male cap (pink)	 \$20
MM-BLU10	LOCKiT extra light male retention cap (blue)	 \$20
MM-GRE10	LOCKiT extended range male retention cap (green)	 \$20
MM-RED10	LOCKiT extra light extended range male retention cap (red)	 \$20
MM-BLA10	LOCKiT male processing cap (black)	 \$20

Green and Red caps are not included in kit.


LOCKiT tools		
MT-LOC10	LOCKiT core tool	 \$79
MT-RSD19	LOCKiT torque wrench driver, short	 \$47
MT-RSD24	LOCKiT torque wrench driver, long	 \$47
MT-LOC11	LOCKiT abutment holder sleeve (5-pack)	 \$45

LOCKiT Abutment kit		
CK-NLC1	LOCKiT abutment kit, H1mm, conical, NP	 \$145
CK-NLC2	LOCKiT abutment kit, H2mm, conical, NP	 \$145
CK-NLC3	LOCKiT abutment kit, H3mm, conical, NP	 \$145
CK-NLC4	LOCKiT abutment kit, H4mm, conical, NP	 \$145
CK-NLC5	LOCKiT abutment kit, H5mm, conical, NP	 \$145
CK-SLC1	LOCKiT abutment kit, H1mm, conical, SP	 \$145
CK-SLC2	LOCKiT abutment kit, H2mm, conical, SP	 \$145
CK-SLC3	LOCKiT abutment kit, H3mm, conical, SP	 \$145
CK-SLC4	LOCKiT abutment kit, H4mm, conical, SP	 \$145
CK-SLC5	LOCKiT abutment kit, H5mm, conical, SP	 \$145
CK-WLC1	LOCKiT abutment kit, H1mm, conical, WP	 \$145
CK-WLC3	LOCKiT abutment kit, H3mm, conical, WP	 \$145
CK-WLC5	LOCKiT abutment kit, H5mm, conical, WP	 \$145
MK-NLC1	LOCKiT abutment kit, H1mm, internal hex, NP	 \$120
MK-NLC2	LOCKiT abutment kit, H2mm, internal hex, NP	 \$120
MK-NLC3	LOCKiT abutment kit, H3mm, internal hex, NP	 \$120
MK-NLC4	LOCKiT abutment kit, H4mm, internal hex, NP	 \$120
MK-NLC5	LOCKiT abutment kit, H5mm, internal hex, NP	 \$120
MK-SLC1	LOCKiT abutment kit, H1mm, internal hex, SP	 \$120
MK-SLC2	LOCKiT abutment kit, H2mm, internal hex, SP	 \$120
MK-SLC3	LOCKiT abutment kit, H3mm, internal hex, SP	 \$120
MK-SLC4	LOCKiT abutment kit, H4mm, internal hex, SP	 \$120
MK-SLC5	LOCKiT abutment kit, H5mm, internal hex, SP	 \$120
MK-WLC1	LOCKiT abutment kit, H1mm, internal hex, WP	 \$120
MK-WLC3	LOCKiT abutment kit, H3mm, internal hex, WP	 \$120
MK-WLC5	LOCKiT abutment kit, H5mm, internal hex, WP	 \$120

 It is possible to use all system components from an existing parallel overdenture solution, like attachments, keys, and caps.

LOCKiT Kit

In order to provide maximum ease-of-use and convenience, the LOCKiT system is packaged as a complete kit that includes everything needed to insure a smooth and convenient procedure. The kit includes the LOCKiT attachment, two silicon rings, metal housing, and plastic caps.

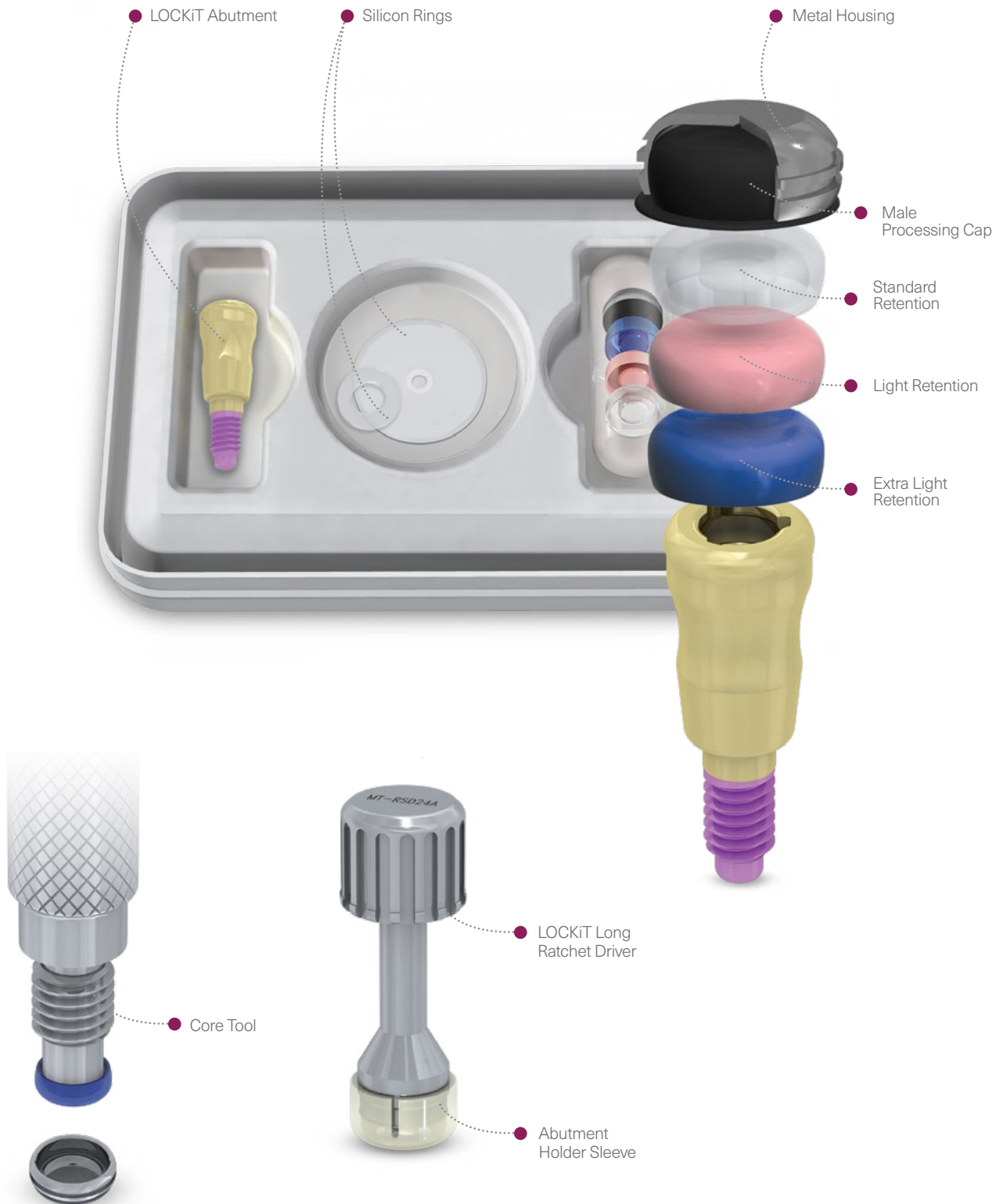



<https://ifu.mis-implants.com>. Adobe Acrobat is required to view the IFU file on the website. This software may be freely downloaded from the Adobe website.

MIS Implants Technologies Inc.
18-00 Fair Lawn Ave. Fair Lawn, NJ,07410
Phone: (201) 797-9144 Fax: (201) 797-9145
Toll Free: 1-866-797-1333
Email: service@misimplants.com
www.misimplants.com

